

Ventilator Cardiac Arrest Procedure and Settings

Paramedic


Cardiac Arrest Settings for Zoll Ventilator

- Choose **Custom** under Start Menu (the following settings will be preprogrammed into the ventilator)
 - Mode: AC (P); this abbreviation indicates Assist Control (Pressure)
 - PIP (Peak Inspiratory Pressure) 25 mmHg
 - BPM (Breaths per Minute) 10
 - FiO₂ 100%
 - PEEP 5
 - Vt (Tidal Volume) Goal around 400 to 600 mL
 - I:E (Inspiratory to Expiratory ratio): 1:2.5



Overview and Device Set-up

CONNECTING THE CIRCUIT AND OXYGEN

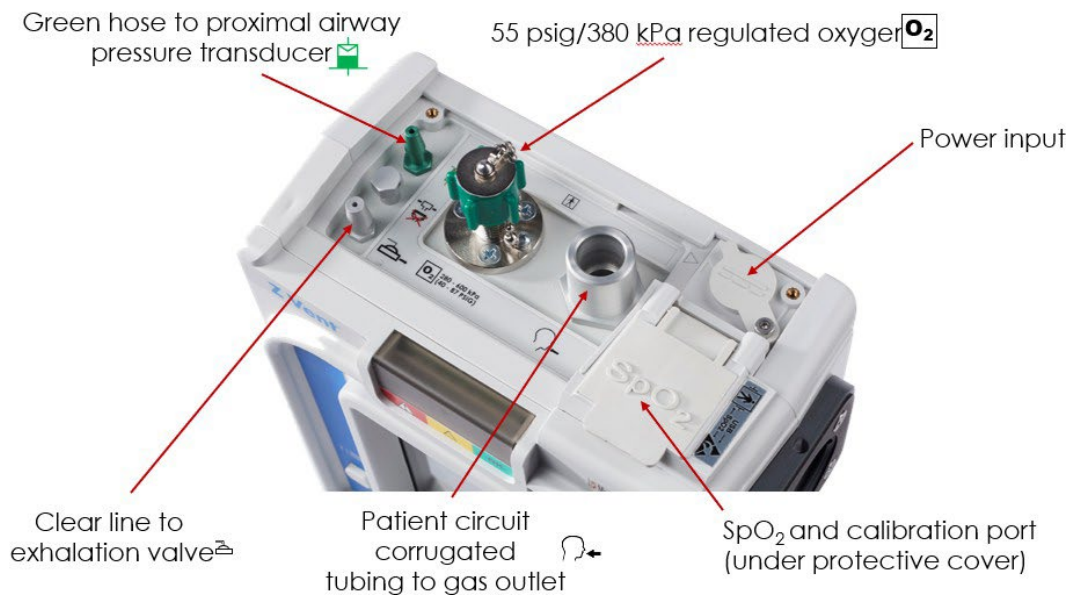

Patient Safety and Protection





Anti-asphyxia protection

The Fresh Gas/Emergency Air Intake acts as the antiasphyxia port, allowing your patient to spontaneously breathe filtered ambient air

Internal Particulate Filter

Part of a 2-stage filtering process.

Add inhalation filter

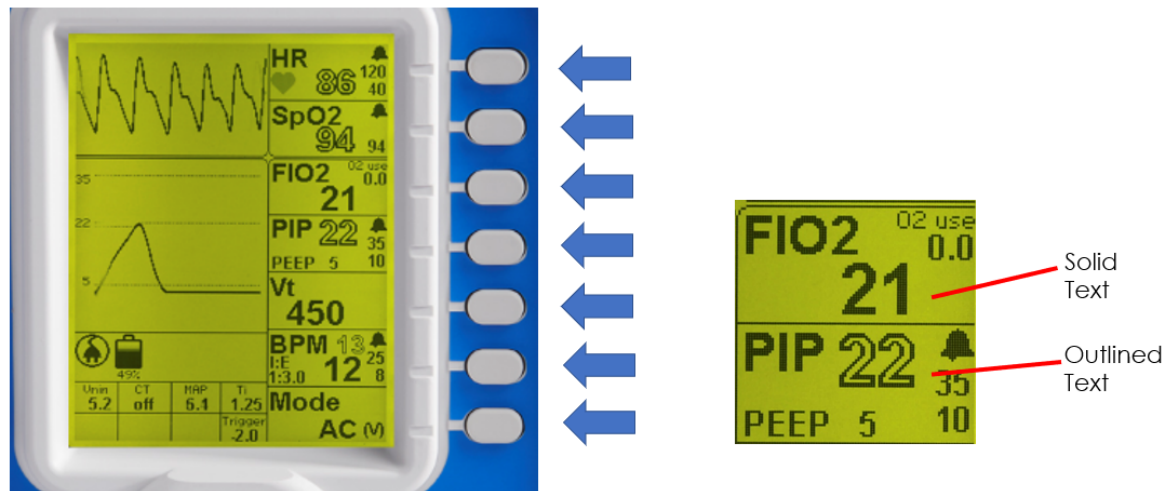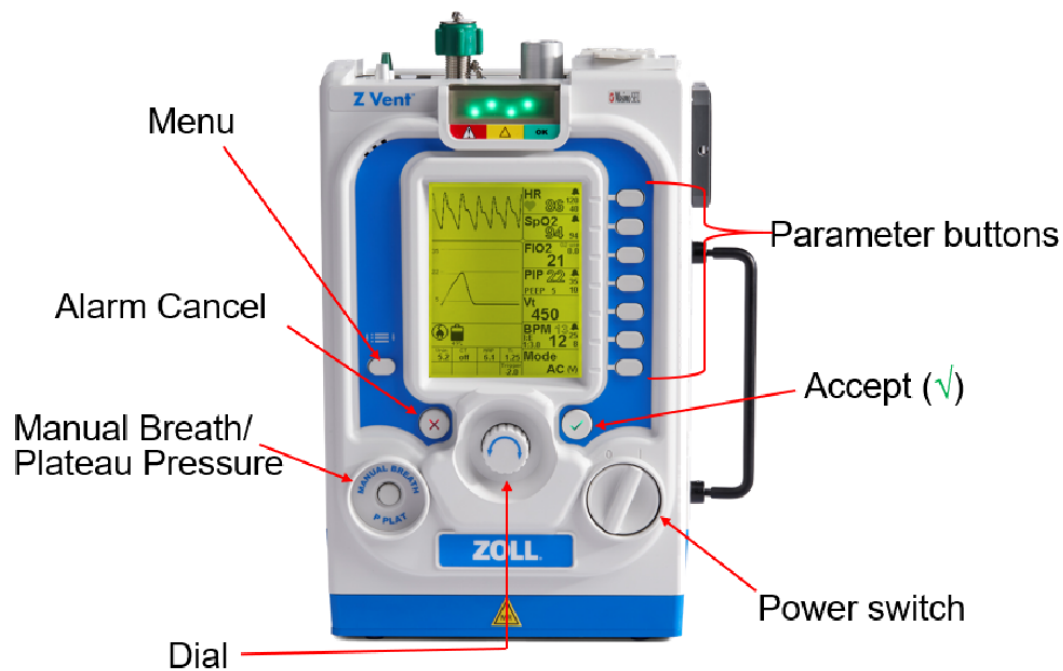

The LCD provides a continuous display of control settings, operating conditions, power, and alarm status information. The display also includes parameter windows with adjacent buttons. These parameters are: Heart rate, SPO₂, FIO₂, Peak Inspiratory Pressure, Tidal Volume, Breaths per minute and Mode of ventilation.

Solid text values can change or be manipulated by the provider. All outlined text values are patient dependent and cannot be changed directly by the provider



- To change a parameter during operation:
- Press button adjacent to parameter window; this highlights the changing parameter.
- Turn Dial to value.
- Press ✓ to confirm the value.



Clinical Use



The Start Menu appears after a quick self-check.

Do not connect patient to the ventilator while the Start Menu is active.

Use the Dial to select a Start Menu option:

Adult Default
Pediatric Default
Mask CPAP Default
Custom
Last Settings

Press ✓ to begin ventilating.

Start Menu – Custom

The Custom option will be used for CPR patients.

Highlight Custom and press Accept (✓)

It is preprogrammed and no setting changes are needed.

Start Menu				HR	120
Adult				stby	40
Pediatric				SpO2	94
Mask CPAP				stby	94
Custom				FIO2	0.0
Last Settings				21	0.0
Auto-PEEP	off			PIP	35
PEEP Leak	on			PEEP	5
PM Due	on			Vt	10
RTC Batt Low	on			450	
Inspiratory Dmnd	on			BPM	25
Intake Restricted	on			12	8
49%				Mode	AC (M)
Umin	0.0	CT	off	MAP	0.0
				Ti	1.25
				Trigger	2.0

Alarm LED Indicators

Green – Patient and ventilator are operating within set parameters.

Yellow – Low Priority: Patient and ventilator are functioning within set parameter.

Red – Medium and High Priority: Patient and/or ventilator are not performing as intended and immediate attention is required.

Pressing Cancel allows a preemptive 30 second alarm mute for low and medium priority alarms.



Alarm Priority

High priority: It is not possible to mute a high-priority alarm. Audio and visual remain until resolved.

Medium priority: It is possible to mute the alarm for 30 seconds only. Audio recurs; visual is constant until resolved.

Low priority: Cancel mutes the audio, visual remains until resolved.

CONNECTING AND DISCONNECTING POWER

- To connect the power supply, locate the power supply connection port on top of the ventilator.
- Open the cover.
- Firmly insert the power cable into the port with triangles aligned. Proper connection will lock the cable into place and prevent accidental removal.
- Give a firm tug to assure the power supply plug has been connected correctly and is locked into place. When the power cord is inserted properly, the user interface will show a power cord in the Icon/Indicator area of the ventilator screen.
- To disconnect, firmly grasp the plug base and pull up. This releases the safety latches, allowing cable removal.
- Never twist or unscrew the power supply cord from the base of the ventilator.